



www.ormesa.com

lenar@ormesa.com

ORMESA



ORMESA®

via A. Da Sangallo, 1 Paciana
06034 Foligno (PG) ITALY
Tel. +39 0742 22927
Fax +39 0742 22637
info@ormesa.com

COMPANY
WITH QUALITY MANAGEMENT
SYSTEM CERTIFIED BY DNV
=ISO 9001:2000=



INGLESE - cod. 22942 - 14-10-09

giuditta

the alternative to being bedridden



THE TILT ADJUSTMENTS of the BACK and (only on the B30 and B12 models) **of the SEAT are continuous and without jerking**

- avoid unbalancing the user
- avoid compression and rubbing the user's column
- Without effort for the assistant

BACKREST and SEAT BOTTOM covered with semi-elastic plastic cord, hard wearing, washable and breathable it avoids the user from sliding forwards. It prevents from bed sores and it is extremely comfortable for those users who are subject to stay for a lengthy time in a stretched or seated position, in alternative to bed

CONTAINED POSTURE without using side supports or belts

GAS SPRING reduces the attendant's effort in raising and lowering and allows a wide range of adjustments

Tyres **PUNCTURE PROOF**


Accessories

828



abdominal belt

834

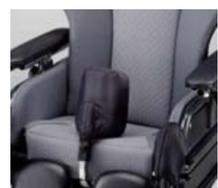

abduction block adjustable in height and backwards-forwards

824-F



fixed table (for Giuditta F30 and F12 models) adjustable backwards-forwards

824-R



reclinable table (for Giuditta B30 e B12 models) adjustable in tilt, backwards-forwards

914



footrest lock Locks the two footrests together to form one unit

892



removable footplate covering

916



padded leg panel in alternative to the calf rests

896



removable tray covering

GIUDITTA F30

FIXED (not tilt-in-space) with 30cm rear wheels



GIUDITTA F12

FIXED (not tilt-in-space) with 4 swiveling 12cm wheels



Comfort & Design

compact **for limited spaces**

easy to handle **for both indoors and outdoors**

comfortable **with low seat for easy transfers**

adjustments **that match physiological needs**

special covering **to limit sliding forward**

resistant **and long-lasting**



PATENTED - Registered Community Design no. 000879424

Contained posture without side supports or belts



holds **those who lack trunk control**

supports **those who tend to slide on one side**

helps **those who seat on their coccyx rather than their buttocks**

accommodates **those who spend a long time lying down or semi-seated**

... why choose



HYGIENIC and BREATHABLE

UPHOLSTERY breathable, fire-retardant, attachable to the seat with buttons easily removable for machine washing

adjustable and



HEADREST multi-directional ADJUSTABILITY

WIDE and wrap-around it is adjustable in height, in a vertical rotating way and forwards-backwards

Giuditta ...



COMFORTABLE


NO-SLIDE COVERING Washable and semi-elastic, it accommodates the user's body shape, maintaining the posture when the backrest is reclined and (only on the B30 and B12 models) when the seat is tilted in any position

highly customizable, adapts to its user



SEAT BOTTOM

depth adjusts in 4 positions (from 40 cm to 45 cm)

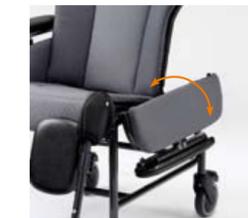

EASY TRANSFERS

LOW SEAT and SWING AWAY SIDE SUPPORTS help the user to get on / off the chair and allow his easy transfers




ARMREST

height-adjustable



NARROW BASE and SWING AWAY LEGRESTS

detachable (only on the Giuditta B30 and B12 models) with FLIP UP FOOTPLATES easy access to elevators and small rooms

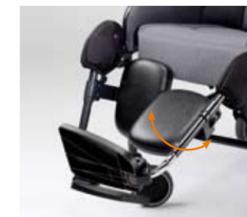

ERGONOMIC LEGREST

adjustable in height and inclination

No-slide, 10° tilt seat



GIUDITTA B30

TILT-IN-SPACE
with 30 cm
rear wheels







GIUDITTA B12

TILT-IN-SPACE
with 4 swiveling
12 cm wheels







Technical Data

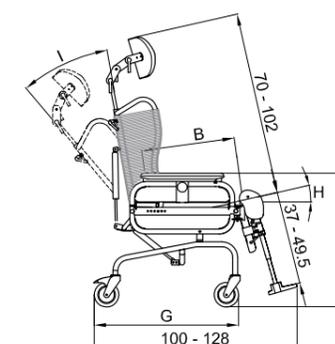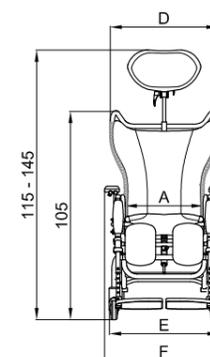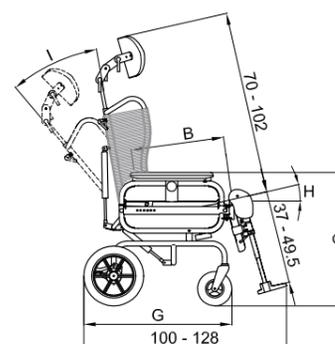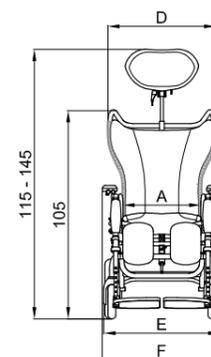

| GIUDITTA B30 (TILT-IN-SPACE) wheels 30 cm and 17 cm | Size 40 | Size 45 | Size 50 |
|--|--|--|--|
| A Seat bottom width | cm 40 | cm 45 | cm 50 |
| B Seat bottom depth | cm 40 - 45 | cm 40 - 45 | cm 40 - 45 |
| C Armrest height | cm 70 - 75 | cm 70 - 75 | cm 70 - 75 |
| D Backrest width | cm 56 | cm 61 | cm 66 |
| E Outside wheel width | cm 59 | cm 64 | cm 69 |
| F Total width | cm 59 | cm 64 | cm 69 |
| G Outside wheel length | cm 76 | cm 76 | cm 76 |
| H Seat bottom inclination (tilt) | 10° - 27° | 10° - 27° | 10° - 27° |
| I Backrest inclination (reclinable) | 35° with B at 40 cm 27° with B at 45 cm | 26° with B at 40 cm 30° with B at 45 cm | 25° with B at 40 cm 30° with B at 45 cm |
| Aid weight | Kg 33 | Kg 34 | Kg 35 |
| Max capacity | Kg 110 | Kg 110 | Kg 110 |

| GIUDITTA B12 (TILT-IN-SPACE) 4 swiveling 12 cm wheels | Size 40 | Size 45 | Size 50 |
|--|--|--|--|
| A Seat bottom width | cm 40 | cm 45 | cm 50 |
| B Seat bottom depth | cm 40 - 45 | cm 40 - 45 | cm 40 - 45 |
| C Armrest height | cm 70 - 75 | cm 70 - 75 | cm 70 - 75 |
| D Backrest width | cm 56 | cm 61 | cm 66 |
| E Outside wheel width | cm 51 | cm 56 | cm 61 |
| F Total width | cm 59 | cm 64 | cm 69 |
| G Outside wheel length | cm 74 | cm 74 | cm 74 |
| H Seat bottom inclination (tilt) | 10° - 27° | 10° - 27° | 10° - 27° |
| I Backrest inclination (reclinable) | 35° with B at 40 cm 27° with B at 45 cm | 26° with B at 40 cm 30° with B at 45 cm | 25° with B at 40 cm 30° with B at 45 cm |
| Aid weight | Kg 30,5 | Kg 31,5 | Kg 32,5 |
| Max capacity | Kg 110 | Kg 110 | Kg 110 |

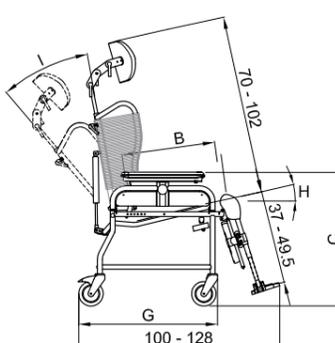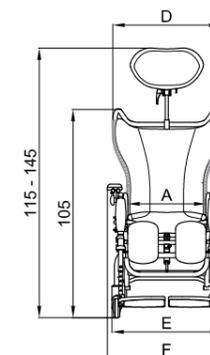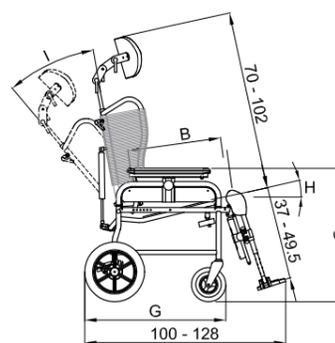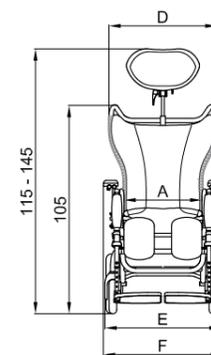

| GIUDITTA F30 (not tilt-in-space) wheels 30 cm and 17 cm | Size 40 | Size 45 | Size 50 |
|--|--|--|--|
| A Seat bottom width | cm 40 | cm 45 | cm 50 |
| B Seat bottom depth | cm 40 - 45 | cm 40 - 45 | cm 40 - 45 |
| C Armrest height | cm 68 - 73 | cm 68 - 73 | cm 68 - 73 |
| D Backrest width | cm 56 | cm 61 | cm 66 |
| E Outside wheel width | cm 59 | cm 64 | cm 69 |
| F Total width | cm 59 | cm 64 | cm 69 |
| G Outside wheel length | cm 70 | cm 70 | cm 70 |
| H Seat bottom inclination | 14° | 14° | 14° |
| I Backrest inclination (reclinable) | 35° with B at 40 cm 27° with B at 45 cm | 26° with B at 40 cm 30° with B at 45 cm | 25° with B at 40 cm 30° with B at 45 cm |
| Aid weight | Kg 29,5 | Kg 30,5 | Kg 31,5 |
| Max capacity | Kg 110 | Kg 110 | Kg 110 |

| GIUDITTA F12 (not tilt-in-space) 4 swiveling 12 cm wheels | Size 40 | Size 45 | Size 50 |
|--|--|--|--|
| A Seat bottom width | cm 40 | cm 45 | cm 50 |
| B Seat bottom depth | cm 40 - 45 | cm 40 - 45 | cm 40 - 45 |
| C Armrest height | cm 68 - 73 | cm 68 - 73 | cm 68 - 73 |
| D Backrest width | cm 56 | cm 61 | cm 66 |
| E Outside wheel width | cm 51 | cm 56 | cm 61 |
| F Total width | cm 59 | cm 64 | cm 69 |
| G Outside wheel length | cm 70 | cm 70 | cm 70 |
| H Seat bottom inclination | 14° | 14° | 14° |
| I Backrest inclination (reclinable) | 35° with B at 40 cm 27° with B at 45 cm | 26° with B at 40 cm 30° with B at 45 cm | 25° with B at 40 cm 30° with B at 45 cm |
| Aid weight | Kg 26,5 | Kg 27,5 | Kg 28,5 |
| Max capacity | Kg 110 | Kg 110 | Kg 110 |