

grillo Accessories



Arm supports **809**
INDOOR & OUTDOOR USE



Weight bars **810**
INDOOR & OUTDOOR USE

The weight bars make the walker heavier according to the size as follows:

- Right + Left size 1 kg. 4
- Right + Left size 2 kg. 4,9
- Right + Left size 3 kg. 5,3
- Right + Left size 4 kg. 5,8



Wheels with adjustable brakes **811**
INDOOR USE ONLY



Directional locks for back wheels **812**
INDOOR USE ONLY



Transparent table **824**
INDOOR USE ONLY



Pair of abduction straps **890**
INDOOR & OUTDOOR USE

GRILLO compensates for lack of balance thanks to its wide base and it allows the User to move by himself thanks to its wheels.

grillo

2007 edition



INGLESSE cod. 22816 - 28/09/09

Grafico: GRAPHIC SERVICE - FOLIGNO

ORMESA®

design and quality made in Italy

Authorized Dealer

CE



ORMESA® S.r.l.

Via A. Da Sangallo, 1 - Loc. PACIANA
06034 FOLIGNO (Pg) ITALY

Tel.: +39 0742 22927 (r.a.)

Fax: +39 0742 22637

E-mail: info@ormesa.com

http://www.ormesa.com

COMPANY
WITH QUALITY MANAGEMENT
SYSTEM CERTIFIED BY DNV
=ISO 9001:2000=

The copyright even if partial is not allowed without the authorization of the manufacturer
The manufacturer reserves the right to make changes without any notice

GRILLO was designed to generate and help the step pattern, to train the pelvis and trunk and to keep their balance, to organize a wider structured vision of space, of distance and of the User's surroundings in his upright position, to develop physical self-sufficiency and psychological independence by experimenting the possibility of moving by oneself.

for **indoor use**



for **outdoor use**

it's available in different sizes and it allows the complete fitness to each user

This Outdoor use model has the same features as the Indoor use, except for the base.



• **CHEST and PELVIC SUPPORTS:** they allow an easy access thanks to the back opening

• **padded PELVIC SUPPORT** for the pelvis. It can be adjusted frontally and backwards to direct the walking with the movement of the barycentre

• **semi-elastic BUFFERS** to protect from any bumps

• **Unpierceable castor WHEELS.** The rear ones have independent brakes. The **accessory 811** on the front wheels helps the adjustment of the flowability of the walking frame, while the **accessory 812** on the back wheels helps its direction

Non-toxic epoxy varnished-plated **FRAME**, with chromium-plated adjustable parts

• adjustable **HANDLEBAR**

it is possible to use only the **PELVIC SUPPORT** if the User has a sufficient trunk control

• **REFERENCES** for the gradual adjustment in height of the chest and pelvic supports

• **LOCKING MECHANISM TO FOLD THE FRAME**

• **padded "ALPINE" PANTIE:** it is comfortable and safe and it holds the User from falling. It is adjustable in height thanks to belts and buckles. It can be easily removed for washing

FOLDING FRAME for transport (except for size 1)

• 15 cm unpierceable castor **WHEELS**

• 25 cm unpierceable fixed **WHEELS**, with **BRAKES**




GRILLO FOR OUTDOOR USE is available in 3 sizes and it allows the use of the walking frame on streets, gardens, etc. without any step and steep slope.



GRILLO for Indoor use

GRILLO 1


Inner circumference of the adjustable pelvic and chest supports, from 47 to 65 cm

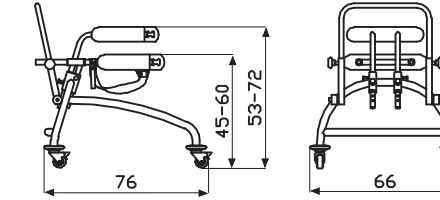

Weight of the aid Kg 13 • Max. load Kg 30

GRILLO 2


Inner circumference of the adjustable pelvic and chest supports, from 56 to 75 cm

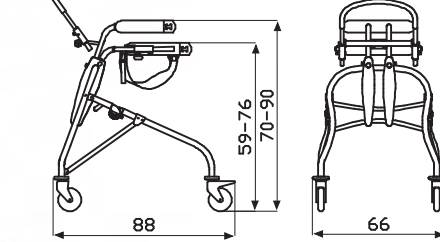

Weight of the aid Kg 18.9 • Max. load Kg 45

GRILLO 3


Inner circumference of the adjustable pelvic and chest supports, from 72 to 92 cm

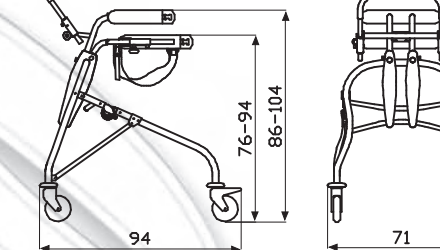

Weight of the aid Kg 20.2 • Max. load Kg 60

GRILLO 4

Inner circumference of the adjustable pelvic and chest supports, from 80 to 100 cm

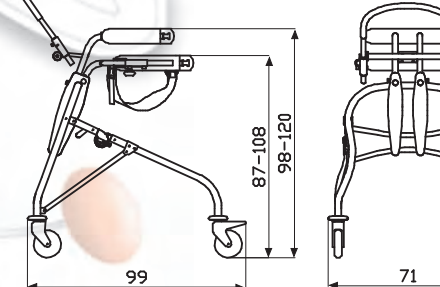

Weight of the aid Kg 20.7 • Max. load Kg 75


Dimensions

in cm and kg

GRILLO for Outdoor use

GRILLO 2


Inner circumference of the adjustable pelvic and chest supports, from 56 to 75 cm

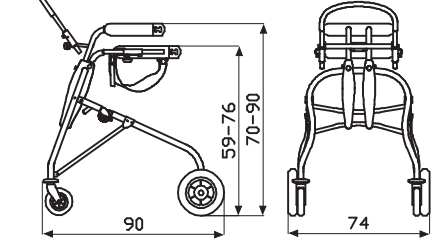

Weight of the aid Kg 20.1 • Max. load Kg 45

GRILLO 3


Inner circumference of the adjustable pelvic and chest supports, from 72 to 92 cm

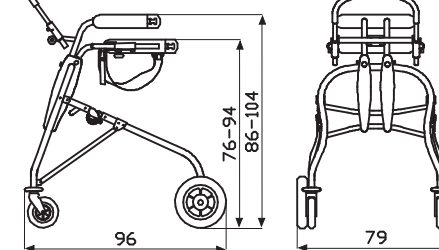

Weight of the aid Kg 21.2 • Max. load Kg 60

GRILLO 4

Inner circumference of the adjustable pelvic and chest supports, from 80 to 100 cm

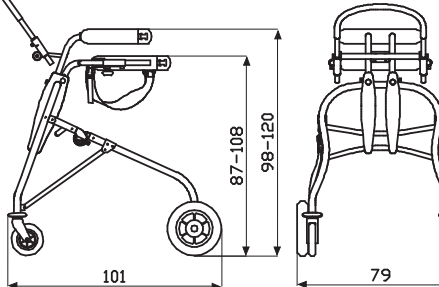

Weight of the aid Kg 21.7 • Max. load Kg 75