


# WASHERS & THERMALATORS

Our Thermalators feature —

- Fully covered bottom promotes easy cleaning
- Thermal protector prevents overheating
- Professional-grade stainless steel construction
- Manufactured from recycled materials 

# Thermalator®



Sturdy stainless steel plastic covered wire rack



Durable handles for easy support

Conveniently-located on/off switch



Locking galvanized steel and rubber casters (on mobile models)



Thermostat allows for precise user-adjustable temperature



Drain valve (on models T-6, T-8 & T-12) speeds up draining for cleanup

Insulated units (T-8 & T-12 models only)



# MOIST HEAT THERAPY

## Thermalator®

The depth of the Thermalator moist heat therapy line offers a unit specific to every kind of treatment setting. Sterile stainless steel stores and maintains Thermal Packs at a consistent heat, which penetrates deep into the muscles to relieve pain and improve circulation.



# NEW!

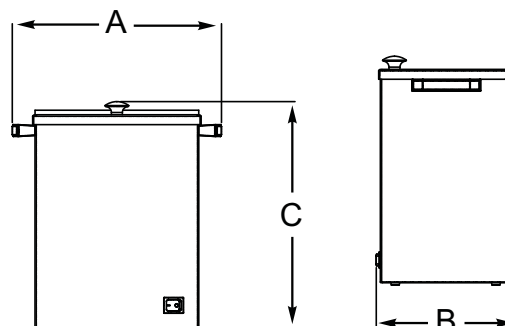
The Little Therm is a miniature of the Thermalator with a four-pack punch. This unit features a clean aesthetic design. Best of all, this product is economically priced, and able to satisfy the most demanding professionals.



Model LT-4  
Little Therm

### Models

heating unit	LT-4
A - length - in (mm)	8-5/8 (219)
B - width - in (mm)	13-1/4 (336)
C - height - in (mm)	14-1/2 (368)
thermal-packs included	4 standard
preview	



- Manufactured from recycled materials.
- Certified per US & Canada Medical Standards; UL 60601-1 & CSA C22.2 No. 601.1

# Thermalator








## Thermalators<sup>®</sup> Moist Heating Units

Whitehall Thermalators maintain a constant water temperature and come standard with exposed heating elements. These units provide for comfortable and effective moist heat therapy, promoting both easy cleaning and safety. The tank body, lid, bottom and tank housing are all made from type 304 stainless steel for years of trouble-free use. The color coded Thermal-Packs fit onto a sturdy stainless steel, vinyl coated wire rack. Inside the unit is a heating element and temperature sensor. Bottom feet keep the unit from slipping and from marring the surface of the countertops. Units have a snap-off thermal protector to prevent overheating. A drain valve located on the back of most units makes them easy to clean.

The Thermalator is offered in sizes that accommodate 4, 6, 8, or 12 hot packs. All units have covered bottoms and are thermostatically controlled. T-4-S, T-6-S and T-8-S are supplied as counter top models. T-8-M and T-12-M are supplied as mobile models with 4 caster wheels (2 are locking).

### Models

heating unit	T-4-S	T-6-S	T-8-S	T-8-M	T-12-M
A - length - in (mm)	10 (254)	15 (381)	16 (406)	16 (406)	16 (406)
B - width - in (mm)	8 (203)	12 (305)	18 (457)	18 (457)	24 (610)
C - height - in (mm)	15 (381)	21 (533)	18 (457)	32 (813)	32 (813)
thermal-packs included	4 standard	1 cervical, 2 oversize, 3 standard	8 standard	8 standard	12 standard
preview					

NOTE: Wire racks, Thermal-Packs and filler hose included with every unit (Filler Hose not included on model T-4-S).

Various combinations and sizes of Thermal-Packs are available.

220V models also available (ET-4-S, ET-6-S, ET-8-S, ET-8-M, ET-12-M).

- Manufactured from recycled materials. 
- Certified per US & Canada Medical Standards; UL 60601-1 & CSA C22.2 No. 601.1 
- Certified per European Standards; IEC 60601-1 



## Moist Heat Therapy Packs

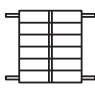

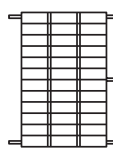
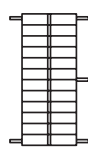
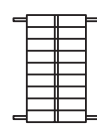

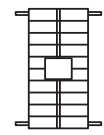

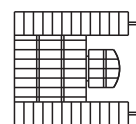
Thermal-Packs contain a filling that evenly distributes and retains heat for up to 30 minutes. Heavy-duty nylon thread keeps them durable; 100 percent cotton makes them comfortable.



### Models TOV1 & TCV1

Oversized & cervical packs

### Models

therapy pack	TST1	TCV1	TOV1	TSL2	TSL1	THF1	TKS1	TMY1	TCB1
length - in (mm)	12 (305)	24 (610)	15 (381)	10 (254)	10 (254)	5 (127)	10 (254)	23 (585)	24 (610)
height - in (mm)	10 (254)	---	24 (610)	24 (610)	18 (458)	12 (305)	20 (508)	---	23 (585)
size	 standard	 cervical	 oversize	 spinal	 spinal	 half	 knee/shoulder	 myofacial	 cervi-back



## Thermal Pack Covers

Terry covers offer even greater comfort and sanitation, and are easily removed from Velcro enclosures to wash and use again.

### Model ATC2



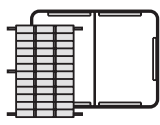
Cervical terry cover with TCV1



### Models ATC4

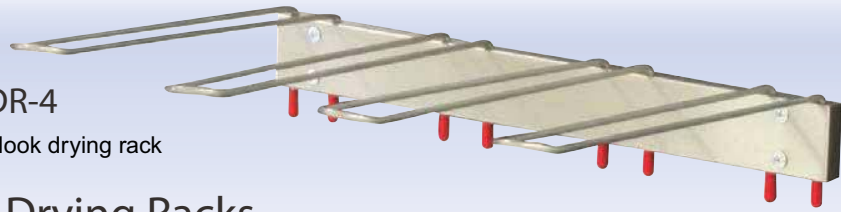
Standard size terry cover

### Terry & Foam Filled units

terry	ATC4	ATC2	ATC6
foam filled	FTC4	FTC2	FTC6
size	standard	cervical	oversize
preview			

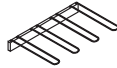
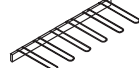


TDR-4  
4-Hook drying rack



## Thermal-Pack<sup>®</sup> Cover Drying Racks

Thermal-pack drying racks are an economical and efficient way to dry Thermal-Pack terry covers. Racks are made from corrosive-resistant stainless steel. Available in four or six hanger configurations.

hooks	4	6
model	TDR-4	TDR-6
preview		

Manufactured from recycled materials



## Thermalator Moist Heating Unit Accessories



### -EG Enviro-Glaze<sup>®</sup>

This coating creates a smooth, protective exterior color finish. Some color choices feature a slight orange peel texture which will not bubble, endures impact and is chip resistant. Enviro-Glaze gives you the durability of stainless steel in a choice of decorative colors.

USAGE: All Thermalators

### -IW4, IW6 Insulating Wrap

This long lasting, commercial-grade grey vinyl cover insulates the Thermalator and cuts energy use in half. It easily installs with Velcro<sup>®</sup> and keeps the exterior surface cool to the touch.

USAGE: IW4 - Thermal Pack T-4-S

IW6 - Thermal Pack T-6-S



### -SBW Side Bottle Warmers

This rack holds (3) 8-ounce dispensing bottles and affixes to the front of the unit to absorb radiated heat in order to warm massage oil or ultrasound lotions.

USAGE: T-4-S & T-6-S Thermalators

Manufactured from recycled materials.



### -STR-1 Side Table Rack

This table is a helpful addition to mobile Thermalators. It comes with three shelves for storing items like bottles or Thermal-Pack covers.

USAGE: T-8-M & T-12-M Thermalators

Manufactured from recycled materials.

